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9. Kontakt

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English

ATTENTION! Read the following instructions before using this product.

1. Product Description

The CAMLOG® Torque wrench is a stainless steel multi-component tool designed to attach to screwdrivers, insertion tools, adapters, and guides. It is supplied with a dual direction wrench mechanism (for insertion/removal). When the selected torque value is reached – a choice of 10, 20, or 30 Ncm – the Torque wrench breaks. This prevents use of an undefined torque value. The Torque wrench can be dismantled into five components for cleaning and care.

2. Indications

Tool for manual thread cutting and manual insertion of CAMLOG® Implants into the bone as far as the definitive position and for the tightening of CAMLOG® Abutment screws and Abutments with a defined torque in Ncm.

3. Contraindications

Any indication not listed here is contraindicated.

4. Safety, Liability, Warranty

The products should be used only by dentists, physicians, and surgeons trained in the implant system. Use of the products requires specialized knowledge and skills in implant dentistry. CAMLOG offers regular courses or technical support for its products. Your local CAMLOG representative will be happy to provide you with information.

Technical support for our products is communicated orally, in writing, with the aid of electronic media and/or by demonstrations. It reflects the state of scientific knowledge and technology current at the time the product was introduced to the market. It does not absolve the user from the obligation to personally test the product for its suitability for the proposed purposes, indications, and procedures. Handling and use of the product are outside our control and are the responsibility of the user. Any liability for damages caused by such actions is excluded.

Use the product only in combination with CAMLOG® Implants and Tools. We do not perform any warranty service or replacement if non-system components have been used. The CAMLOG® Implant System is part of an overall concept and is to be used only with original system parts and tools in accordance with the recommendations given in the instruction manuals and the directions of the manufacturer. All components of the CAMLOG® Implant System are optimally matched to one another.

The use of non-system components may interfere with system function and will lead to the exclusion of all guarantee and replacement services by ALTATEC GmbH/CAMLOG.

Tools and system components are designed for specific implant and prosthetic lines and implant diameters. Their use with different implant or prosthetic lines or different implant diameters can lead to mechanical malfunction of system components, tissue damage, or esthetically unsatisfactory results. For this reason, some implant and prosthetic lines are given to some extent their own components with matching tools. Pay attention to the color-markings when choosing tools for the required implant diameter.

As part of our sales and delivery terms, we guarantee the flawless quality of our products. Information on existing patents, brands, or other intellectual property laws are not legally binding.

As a general rule, all products must be protected from the risk of aspiration and swallowing during intraoral use. With regard to the indications listed for the particular products, take note that any indication that is not listed is contraindicated.

Unless they are specifically indicated as sterile, products of the CAMLOG® Implant System are delivered non-sterile.

Our products are sold only to physicians or dentists or on their prescription. Not all parts are available in all countries.

5. Use

5.1 Attachment of Prosthetic Components


For definitive installation in the implant, all abutment and prosthetic screws as well as abutments must be tightened to the preset torque with the Torque wrench.

WARNING!
To obtain maximum preload values for the screws, abutments must be retightened after approx. 5 minutes with the same torque. This will help prevent screw loosening.

Tool Installation


The tools lock in during their insertion in the wrench wheel.

- For the insertion function: “IN” is visible on the wrench button
- For the removal function: “OUT” is visible on the wrench button.
- To remove, press on the tool button with the finger and remove downward. Do not use pliers!

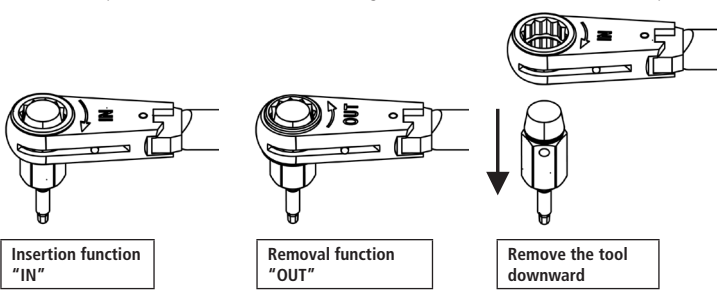
Setting Torque Limits

The release torque (10, 20 and 30 Newton centimeters) is set by screwing in or screwing back the torque setscrew.

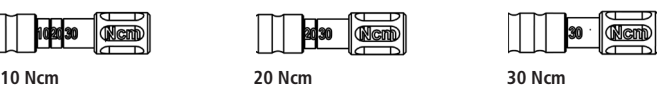The Torque wrench breaks when the release torque is reached. After that, do not tighten further!

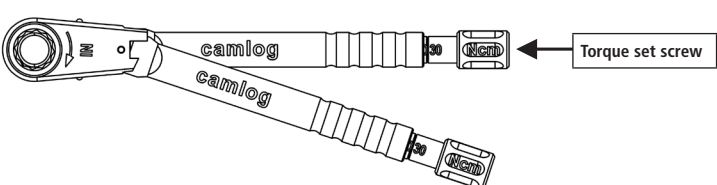Torque settings for the Torque wrench	
Implant cover screw	Hand tightened
Gingiva former	Hand tightened
Impression post	Hand tightened

Abutment and Prosthetic Screws:

– Provisional abutment	Hand tightened
– Ceramic abutment	20 Ncm
– Esthomic® abutment	20 Ncm
– Telescopic abutment	20 Ncm
– Universal abutment	20 Ncm
– Gold-plastic abutment	20 Ncm
– Logfit™ abutment	20 Ncm
– Prosthetic screw for bar abutment	15 Ncm

Abutments:	Ø 3.3 mm	Ø 3.8 – 6.0 mm
– Bar abutment	20 Ncm	30 Ncm
– Ball abutment	20 Ncm	30 Ncm
– Locator® abutment	20 Ncm	30 Ncm

NOTE: In the laboratory, abutments should not be tightened with the Torque wrench! The listed values are relevant only for the clinical setting.

WARNING! To ensure optimal preloading, NEW abutment screws should be used for the definitive attachment of the abutment.

5.2 Manual Thread Cutting and Implant Insertion

For this procedure the wrench is used in its locked position. The lock is set as follows:

- Loosen the torque setscrew (counterclockwise).


- Slide back the handle.


- Rotate handle 90°.

- Slide handle forward to the locking position in the wrench head.

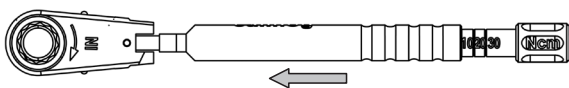- Fasten the handle with the torque setscrew (clockwise). The torque function is now locked.

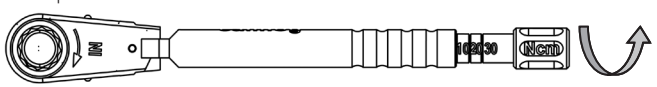
Torque function locked

6. Cleaning, Sterilization, Maintenance

6.1 Disassembly


The Torque wrench should be disassembled after use into its five components.

- To do this, completely unscrew the torque setscrew (counterclockwise) and pull out the spring and handle.
- Pull back the fastening pin with your thumb and forefinger and remove the wrench wheel from the wrench head.
- Clean the individual parts with a small plastic brush under running water and place them in a disinfectant-cleaning liquid (for the duration, see the manufacturer's instructions).
- Dry all parts.
- Lubricate the indicated places () with contra-angle handpiece oil or spray and reassemble the Torque wrench.

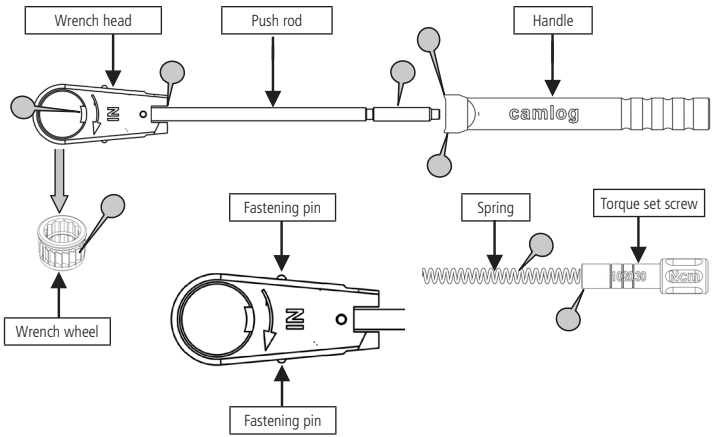6.2 Assembly

- First install the wrench wheel: with your thumb and forefinger, pull back the fastening pin on both sides and insert the wrench wheel narrow side in the wrench head (“IN” must show).
- Slide the handle back on, insert the spring in the handle, insert and tighten the torque setscrew (clockwise).
- Sterilize the Torque wrench in the autoclave prior to clinical use (see “Preparation Instructions for the CAMLOG® Implant System”, Art. No. J8000.0032).

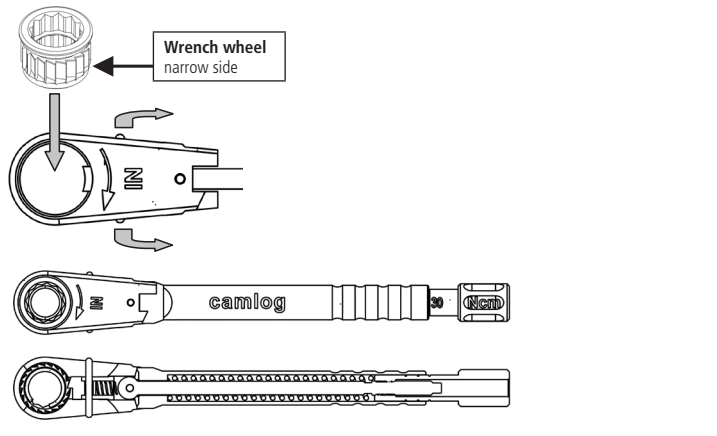WARNING! Non-sterile shipped CAMLOG® products must not be sterilized in the CAMLOG original packaging!

7. Further Documentation

Further information on our products is available at www.camlog.com.

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